

Dermatology

Derma Pen



Derma Pen Basic with 0.25mm to 2.0mm



Derma Pen Pro with 0.25mm to 2.0mm




Derma Pen Ultima with 0.25mm to 2.5mm

Salient Features

The Dermapen is a motorized surgical instrument that uses a disposable needle cartridge with 9 / 12 / 36 micro-needles.

- No cross infection / One-time disposable needle consumables
- Speed – faster than manual applications
- Short session times
- Adjustable needle depth from 0.25 to 2 / 2.5 mm depending on the model
- Minimal pain
- Lightweight handpiece
- Easy to use
- Low cost consumable
- User independent

Dermapen Micro-Needling Depth Chart

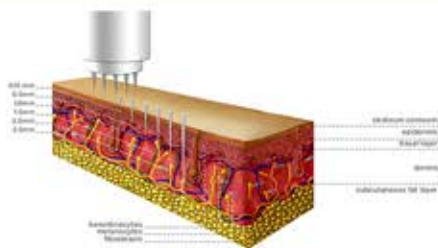

REGION	THIN SKIN	THICK SKIN
Forehead	0.25-0.5mm	0.5-0.75mm
Between Eyebrow	0.25-0.5mm	0.5-1.0mm
Nose	0.25mm	0.5mm
Eye Area	0.25mm	0.25-0.5mm
Cheek Bone	0.5mm	0.5-1.0mm
Cheek and Chin	0.5-1.0mm	1.0-2.0mm
Upper Lip	0.25mm	0.25-0.75mm
Facial Scarring	1.0-1.25mm	1.5-2.0mm
Stretch Mark	1.5-2.0mm	1.5-2.0mm
Surgical Scar Removal	2.0-2.5mm	

Reference Chart



Functions

- Anti-Puffiness
- Stretch Mark Removal and Stretch Mark Treatment
- Acne Scar Removal, Acne Scar Treatment
- Anti Wrinkle Treatment, Anti Aging
- Hair Loss Treatment or Hair Restoration
- Surgical Scar Removal.